



**IMPLANTS** DIFFUSION INTERNATIONAL



**ID<sup>CAM</sup>ST**

**FABRIQUÉ EN**  
**FRANCE**

**French** developer and manufacturer since 1987

# ID<sup>CAM</sup> ST / STANDARD

## Range

### PRESENTATION

The ID<sup>CAM</sup> ST range type IM benefits from the SLA + TiO<sub>2</sub> state of surface initiated by IDI and used since 1987.

This self-condensing implant stands 75 N.cm screwing stress without being affected. The ID<sup>CAM</sup> ST implant draws special attention to itself due to its Switching Cone neck and to its cylindro-tapered body identical to a dental root.

The angulation, the state and the depth of the threads are specially studied to optimize the primary stabilization in any bone density and favour the immediate loading.

### Implant ID<sup>CAM</sup> ST features:

- Cylindro-tapered-shaped
- 2.5° morse taper
- Titan alloy TAL6VELI
- SMA + TiO<sub>2</sub> state of surface
- Cam retention
- Switching Cone
- Anti-unscrewing grooves
- Progressive and condensing threads
- Convex apex

Switching  
Cone

Condensing  
high threads

Anti unscrewing  
groove

Self-drilling  
low threads

CONVEX  
APEX


2,5° Morse Taper

Cam connection

Internal Universal  
thread ISO

Morse taper  
and identical  
connection for  
all ID<sup>CAM</sup> range.

# MORSE TAPER CONNECTION & TUBE'N'TUBE

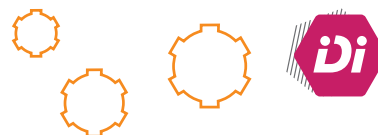

## IDCAM ST IMPLANTS RANGE

| DESCRIPTION | LENGTH<br>Color code* | DIAMETER      | PLATFORM | P/N       |
|-------------|-----------------------|---------------|----------|-----------|
| STANDARD    | 8 mm ●                | <b>3,7 mm</b> | 3,6 mm   | IDCS0835  |
| STANDARD    | 8 mm ●                | <b>4,2 mm</b> | 3,6 mm   | IDCS0842  |
| STANDARD    | 8 mm ●                | <b>5,2 mm</b> | 3,6 mm   | IDCS0852  |
| STANDARD    | 10 mm ●               | <b>3,7 mm</b> | 3,6 mm   | IDCST1035 |
| STANDARD    | 10 mm ●               | <b>4,2 mm</b> | 3,6 mm   | IDCST1042 |
| STANDARD    | 10 mm ●               | <b>5,2 mm</b> | 3,6 mm   | IDCST1052 |
| STANDARD    | 12 mm ●               | <b>3,7 mm</b> | 3,6 mm   | IDCST1235 |
| STANDARD    | 12 mm ●               | <b>4,2 mm</b> | 3,6 mm   | IDCST1242 |
| STANDARD    | 12 mm ●               | <b>5,2 mm</b> | 3,6 mm   | IDCST1252 |
| STANDARD    | 15 mm ●               | <b>3,7 mm</b> | 3,6 mm   | IDCST1535 |
| STANDARD    | 15 mm ●               | <b>4,2 mm</b> | 3,6 mm   | IDCST1542 |
| STANDARD    | 15 mm ●               | <b>5,2 mm</b> | 3,6 mm   | IDCST1552 |

\*On each implant packaging there is a small colored sticker to match with the implant length. The code for each color is related to the one found on the RBS conical drills and TURBOdrill® for IDCAM:

- Length 8 mm
- Length 10 mm
- Length 12 mm
- Length 15 mm

All the dimensions are in millimeters.

## IMPORTANT CONSIDERATIONS:




The 3,5mm diameter implants must be used only for lower and upper lateral incisors.

**The IDCAM implants are supplied with a closing and a healing cap.**


**To improve aesthetic quality, it is recommended to position the implant 1mm under the bone crest.**

It is recommended to set as many implants as lacking natural roots in the patient mouth in order to secure the lasting of the prosthesis.


The length choice and diameter of implant must be based on the density determined by CT-scan.

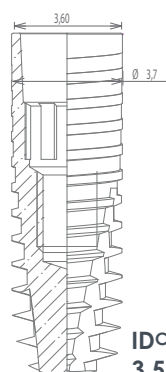

IDCAM ST  
3,5

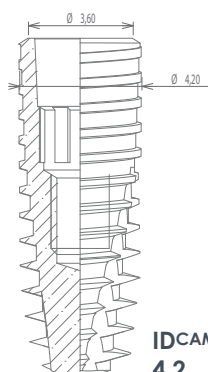

IDCAM ST  
4,2

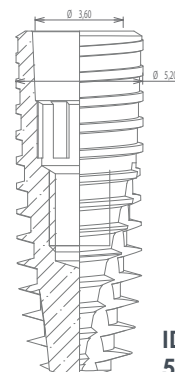

IDCAM ST  
5,2

# ID<sup>CAM</sup>ST SURGICAL PROTOCOL

## DRILLING PROTOCOL WITH A ID<sup>CAM</sup> TURBODRILL SURGICAL SET (PROTOCOL WITHOUT COLLECTING BONE)

Example for an Ø4,2 mm implant

### • IN THE MAXILLA

1. Use the pilot drill of Ø2 mm at 650 rpm with irrigation;
2. Use the 2nd drill of Ø2,8 mm between 1000 and 1500rpm with ample irrigation;
3. Use the 3rd drill of Ø3,5 mm between 1000 and 1500rpm with ample irrigation;
4. Use the 4th drill of Ø4,2 mm between 1000 and 1500rpm with ample irrigation;  
- Four-blade drills must be used without back-and-forth movement.
5. Then screw the implant 1 mm under crestal.

### • IN THE MANDIBLE

1. Use the pilot drill of Ø2 mm at 650 rpm with irrigation;
2. Use the 2nd drill of Ø2,8 mm between 1000 and 1500rpm with ample irrigation;
3. Use the 3rd drill of Ø3,5 mm between 1000 and 1500rpm with ample irrigation;
4. Use the 4th drill of Ø4,2 mm between 1000 and 1500rpm with ample irrigation;
5. For **dense bone**, pass the C1840 drill, with depth markings 8,10,12,15, and 18 mm, to the corresponding implant length, between 1000 to 1500 rpm under ample irrigation.
6. Then screw the implant 1 mm under crestal.

**Expert advice:** After the last drilling, rinse the socket abundantly with physiological serum before implant placement.

### • IMPLANT PLACEMENT

1. Insert the implant at 1.0 mm under-crestal using a contra-angle with dental shank screwdriver P/N: 1046 or 1146, at 35 rpm with a torque of 30 to 40 N.cm, or manually using the ratchet P/N: 415 and the screwdriver P/N: 0046 or 0146.
2. Screw the closing cap and the healing cap manually at 5 N.cm

**Expert advice:** Soak the implant, the closing cap and the healing cap in an antibiotic (tobramycin or gentamicin type – 75 mg) diluted in 20 cl of physiological serum before placing it.

**Life cycle :** The drills shall be replaced after 15 to 20 uses or as soon as their cutting capacity decreases. Used drills must be decontaminated or treated as **DASRI (Infectious Risk Care Activity Waste)**.

### **Drilling with bone recovery (without irrigation):**

Use the drill sequence of progressive diameters at 100 rpm without irrigation.

Soak the drill, loaded with bone, in a capsule filled with physiological serum, so that the bone detaches and settles at the bottom. Aspirate serum and recover bone ready for bone graft.



IMPLANTS DIFFUSION INTERNATIONAL

23/25 rue Émile Zola - 93100 Montreuil - FRANCE

Tél. : +33 (0)1 48 70 70 48 – Email : info@idi-dental.com

**Find all IDI implant ranges on our website:**  
**[www.idi-dental.com](http://www.idi-dental.com)**