

# The ID'BOX


## Supporting your surgery

### ID'PARALLELISOR

FOR MORE PRECISION, THIS ID'BOX GUIDES YOUR IMPLANT SETTING.

#### PRESENTATION

This box is indicated for the placement of the implants in a paralleling way to maximize the life of the prosthetic attachments.

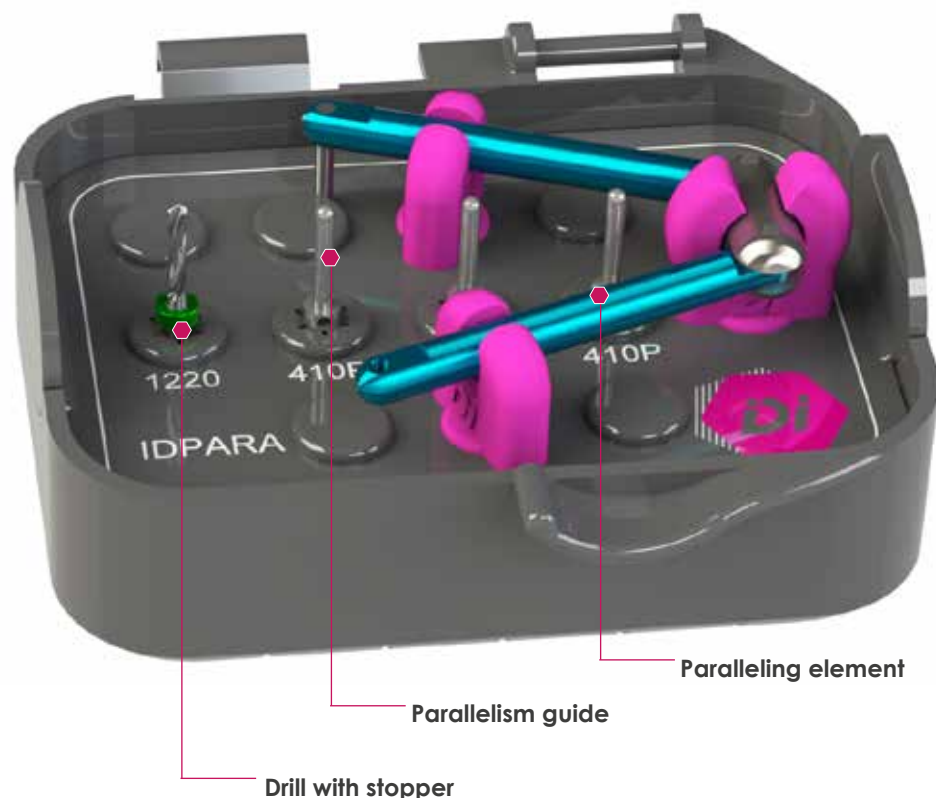

#### CONTENTS :

- 1 drill with stopper
- 3 parallelism guides
- 1 paralleling element (articulated instrument)

#### ID-PARALLELISOR PROTOCOL

##### STEP 1

Drill perpendicular to the basal bone with the drill (green stopper).



##### STEP 2

Insert the stem of the paralleling element into the freshly drilled socket. Shift the parallelisor axis and drill through the hole located in the parallelisor arm where the 2nd implant must emerge.

